

*No. 1 choice
of hospitals
& mothers*

Optimise hospital processes with disposable products




- I No need for cleaning prior to first use
- I Single- or one-day-use
- I Minimisation of breastmilk contamination risk
- I Convenient and simple: Disposable after use



Improve your hospital process – with our simple and disposable products

Our breastfeeding products are the result of more than 30 years of evidence-based research. Thanks to this accumulated knowledge, Medela is the partner for complete solutions in the hospital sector. Our pump sets, breastmilk bottles and feeding solutions help to get precious breastmilk safely to where it is needed: in the infant's tummy. After all, we all want the same thing: happy and healthy mothers and infants.

With its disposable product range, Medela greatly simplifies the daily routine involved in breastfeeding. All these products can be used straight away without cleaning prior to first use. This supports hygiene and a safe use. Cost savings are likely, due to the elimination of the time, cost and energy-intensive sterilisation/disinfection process.

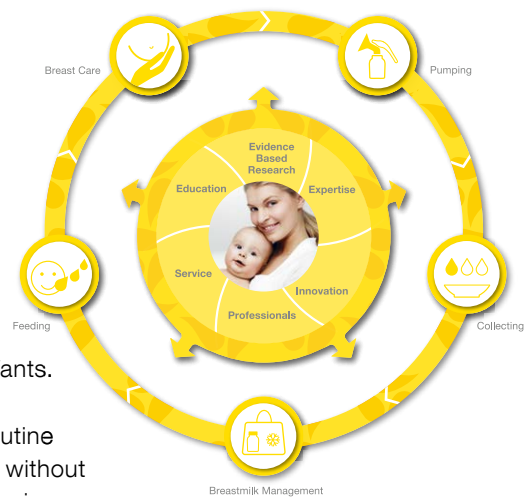

Breastfeeding Solution Circle

Medela pump sets: Hygienic and convenient handling

Our pump sets are designed for optimised and convenient handling. Every mother gets a new personal set to minimise the risk of breastmilk contamination. Since the One-Day Pump Sets can be used for up to eight pumping sessions within 24 hours, one set a day is enough. During this time the mother takes responsibility for her pump set (e.g. cleaning, storing). The simplicity of the design makes it easy for her to clean and handle. This saves you costs and time which you can devote to your other duties and attending to other patients. To satisfy the mothers' individual needs we offer different breastshield sizes to maximise comfort and efficiency while pumping. Moreover, the appropriate bottle size can be chosen according to the personal milk output.

Medela breastmilk bottles: Right size for every use

Medela Disposable Bottles are used to collect, store, portion and feed breastmilk. In a milk bank, the bottles can also be used for pasteurisation. We provide different sizes of bottles: 35ml, 80ml, 150ml and 250ml.

The 35ml Colostrum Container is a great example of how all Medela Disposable Breastmilk Bottles are designed with healthcare professionals and mothers in mind. The curved bottom minimises the loss of colostrum or breastmilk during transfer to a syringe. The smaller size helps to keep mothers motivated to pump by setting realistic expectations of initial milk production.


The breastmilk bottles are graduated which enables easy checking of the content. In addition, you can note the infant's name, the volume expressed, as well as the time and date of expression on the additionally available label.



One-Day Pump Sets



Disposable Breastmilk Bottles

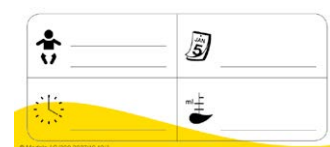

Label for Disposable Breastmilk Bottles

Natural feeding experience: Medela differentiated feeding solution

Calmita is a research-based hospital feeding solution that allows newborns to train and apply their natural sucking behaviour. The integrated vacuum-controlled valve allows the infant to decide when to drink and when to pause. Calmita is available in two versions in which the threshold levels of the valve are different. Milk flows when the infant reaches a certain vacuum. The neonate creates its own sucking rhythm, thereby efficiently removing just the right amount of milk at an individually suitable pace.



Calmita Starter



Calmita Advanced



Differentiated disposable feeding solution

Conventional way: Medela conventional teats

Our individually packed conventional disposable teats help to simplify feeding an infant with breastmilk. Two different milk flow rates accommodate the different needs of newborn infants. As an alternative to teat feeding, Medela offers a wide range of special feeding devices, for example devices to support infants with cleft lip/palate abnormalities.



Standard



Premature

Conventional disposable teats



Medela breastmilk bottles and all our products that come in contact with breastmilk are and always have been BPA-free.



Medela Ready-to-Use products are hygienically safe to use without prior cleaning



Medela sterile products are validated according to the applicable standards for sterile medical devices

Medela's reusable products

For hospitals with traditional cleaning and sterilising practices.



Pump sets

- I Designed for use by multiple mothers
- I Suitable for autoclaving
- I Single and double pump sets available
- I Various PersonalFit breastshield sizes available (S, M, L, XL, XXL)

Breastmilk bottles

- I Designed for use by multiple mothers
- I Suitable for autoclaving
- I Sizes 80ml and 150ml available

Standard teats

- I Small flow for preterm infants
- I Medium flow for term infants
- I Suitable for autoclaving



Symphony Pump Set



Basic Kit



80ml



150ml



Slow flow



Medium flow

www.medela.com

Disposable products

Pump sets

Article number	Product	Materials ¹⁾	Qty
008.0176	Symphony One-Day Pump Set with PersonalFit Size M (24mm)	PP, TPE	54
008.0177	Symphony One-Day Pump Set with PersonalFit Size L (27mm)	PP, TPE	45
008.0178	Symphony One-Day Pump Set with PersonalFit Size XL (30mm)	PP, TPE	45
008.0179	Symphony and Lactina One-Day Pump Set with PersonalFit Size M (24mm)	PP, TPE	72
008.0180	Symphony and Lactina One-Day Pump Set with PersonalFit Size L (27mm)	PP, TPE	60
008.0181	Symphony and Lactina One-Day Pump Set with PersonalFit Size XL (30mm)	PP, TPE	60

Breastmilk bottles

200.2110	35ml Colostrum Container	PP, PBT	100
200.2957	80ml Disposable Breastmilk Bottles	PP, HDPE	480
200.2956	150ml Disposable Breastmilk Bottles	PP, HDPE	320
200.4802	250ml Disposable Breastmilk Bottles	PP, HDPE	160

Teats

008.0257	Calmita Starter	PP, TPE	260
008.0258	Calmita Advanced	PP, TPE	260
200.3200	Disposable Hospital Teat, Standard, Medium Flow	PP, TPE	100
200.3201	Disposable Hospital Teat, Premature, Slow Flow	PP, TPE	100

Reusable products

Pump sets

800.0557	Symphony Hospital Multibox with PersonalFit Size M (24mm)	PP, Silicone	20
200.1141	Symphony Hospital Multibox with PersonalFit Size L (27mm)	PP, Silicone	20
200.1143	Symphony Hospital Multibox with PersonalFit Size XL (30mm)	PP, Silicone	20
800.0605	Basic Kit (Breastshield M (24mm), valve, bottle, lid)	PP, Silicone	25

Pump set accessories

200.3887	Label for Disposable Breastmilk Bottles	-	5x2'500
008.0182	Label Dispenser	-	1
800.0608	Symphony Membrane Cap and Protective Membrane	-	25
800.0609	Lactina Piston/Cylinder incl. Cylinder Holder	-	5
800.0659	Silicone Tubing for Symphony and Lactina	-	25
200.0571	PersonalFit Breastshield Size S (21mm)	PP	25
800.0668	PersonalFit Breastshield Size M (24mm)	PP	25
800.0711	PersonalFit Breastshield Size L (27mm)	PP	25
800.0712	PersonalFit Breastshield Size XL (30mm)	PP	25
800.0854	PersonalFit Breastshield Size XXL (36mm)	PP	25

Breastmilk bottles

810.0137	80ml Reusable Breastmilk Bottle	PP	1
800.0625	150ml Reusable Breastmilk Bottle	PP	25
800.0627	Lid	PP	25
800.0628	Disc	PP	25

Teats

008.0119	Reusable Hospital Teat, Slow Flow	Silicone	25
008.0120	Reusable Hospital Teat, Medium Flow	Silicone	25

Medela also offers Hospital Grade Breastpumps: Symphony and Lactina. Find out more at www.medela.com or contact your local Medela representative.

¹⁾ in contact with breastmilk

